

2015 EDITION Ⅲ



# SAIKANG MEDICAL

# JIANGSU SAIKANG MEDICAL EQUIPMENT CO.,LTD. JIANGSU SAIKANG MEDICAL CO.,LTD.

Addr.: No. 35 Lehong Road, Modern Agriculture Demonstration Park, Zhangjiagang City, Jiangsu Province, China.

Tel.: 0086-512-58525850

Fax.: 0086-512-58520713 Mob.: 0086-13962479538

P.C.: 215600 Web.: www.chinasaikang.com www.jssaikang.com

E-mail: inquiry@jssaikang.com www.jssaikang





Company profile

convenient transportation access.

We are combined with manufacture, service, and supply of Hospital Furniture, important care on our customer's satisfaction. So far, over 80 countries and philosophy of being firmly committed to "Integrity is the base. Customer is the first".

Managed by a professional team who has a wide range of managerial experience. great enthusiasm and efficacy. Our vision we provide the world-class products in the near future.

Certificate







Emergency, Surgical products. We are devoted to providing a one-stop sales model to customer; We give more regions have enjoyed our products and service. Our success has been mainly achieved by following the corporate

our Sales-team, Research and Development department, Production department, QC department are working with and service to improve health-care system worldwide to the top level. Wish you can enjoy our products and service











### CONTENTS

Ambulance Stretcher Trolley 01 Stretcher Base & Accessories 05

Stretcher 07

Stretcher Accessories 17

Special Stretcher 19 Soft Stretcher 22

Resuscitator / Face Mask 24

Oxygen Intake Device 25

Defibrillator & Tent Bed 26

Emergency Training Model 27

Mortuary Refrigerator 28

Body Bag / Church Trolley 29

Jiangsu Saikang Medical Equipment Co.,Ltd. has been established since 2002. It is located in Zhangjiagang City, along the Economic circle of Yangtse River Delta, enjoying the highly developed economic expansion and convenient transportation

Our success has been mainly achieved by following the corporate philosophy of being firmly committed to "Integrity is the base, Customer is the first". Our vision is "Race against time, healthier your life.", to provide the world-class products and services to improve health-care system worldwide to the top level. Our mission is to achieve 60 percent of the market share in our products range.



# Stretcher Trolley For Ambulance-







Wodel	(High Position)	(Low Position)	Max. Aligie of the back	Load Bealing	14.44.
SKB039(A)	1900*550*500mm	/	75°	159kg	25.5kg
SKB039(B)	1900*550*500mm	1	75°	159kg	25.5kg
SKB039(D)		1900*550*250mm	75°		













Model	Product Size (High Position)	Product Size (Low Position)	Max. Angle of the back	Load Bearing	N.W.
SKB039(C)					
SKB039(E)	1960*550*860mm	1960*550*250mm	85°	159kg	34kg





# Stretcher Trolley For Ambulance-



Model	(High Position)	(Low Position)	Max. Angle of the back	Load Bearing	N.W.
SKB039(G)			75°		





Model	Product Size (High Position)	Product Size (Low Position)	Max. Angle of the back	Load Bearing	N.W.
SKB040(A)	1900*550*880mm	1900*550*540mm	1	159kg	34kg
		1900*550*54 0mm	65°		



### Stretcher Base Accessories





SKB3A204-1 SKC-3D4 Stretcher Base





I.V. Stand





### Foldable Stretcher



## SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2140\*540\*120mm Folded Size(L\*W\*H): 580\*210\*140mm



## SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2080\*550\*130mm Folded Size(L\*W\*H): 1040\*170\*90mm



Folded Size(L\*W\*H): 610\*300\*120mm



# SKB1A01-2

SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2150\*550\*150mm Folded Size(L\*W\*H): 550\*260\*170mm



SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2150\*550\*130mm Folded Size(L\*W\*H): 1820\*90\*170mm





SKB1A03

SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 1850\*500\*180mm Folded Size(L\*W\*H): 960\*540\*160mm

SKC-3A9 Foldable Stretche

Full Size(L\*W\*H): 2080\*550\*130mm Folded Size(L\*W\*H): 1040\*170\*90mm



### SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 1850\*500\*240mm





### Foldable Stretcher



# SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2300\*170\*550mm Folded Size(L\*W\*H): 1030\*130\*230mm



# KC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 1850\*500\*200mm Folded Size(L\*W\*H): 925\*50\*120mm



### SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 1850\*500\*30mm Folded Size(L\*W\*H): 925\*500\*70mm



SKC-3A9 Foldable Stretcher With Body Ba

Full Size(L\*W\*H): 2080\*530\*350mm Folded Size(L\*W\*H): 1080\*530\*120mm







### SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 1850\*500\*30mm Folded Size(L\*W\*H): 930\*500\*70mm



### KC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2200\*550\*150mm Folded Size(L\*W\*H): 950\*200\*90mm





### SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2200\*550\*150mm Folded Size(L\*W\*H): 950\*200\*100mm



### Foldable Stretcher Stair Stretcher







SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 1880\*450\*40mm Folded Size(L\*W\*H): 480\*450\*150mm



# SKC-3A9 Foldable Stretcher

Full Size(L\*W\*H): 2200\*550\*150mm Folded Size(L\*W\*H): 1800\*120\*150mm



SKC-3A9 Foldable Stretche









SKB1C01-1	900*520*1360mm	1120*520*200mm	159kg	
			159kg	



# SKC-2B3 Stair Stretcher

Full Size: 690\*510\*900mm Folded Size: 180\*510\*900mm



#### SKB1C04 SKC-2B3 Stair Stretcher Full Size: 690\*510\*900mm

Folded Size: 180\*170\*590mm



### SKB1C05

SKB1C02-1

SKC-2B3 Stair Stretcher

SKC-2B3 Stair Stretcher Full Size: 690\*510\*900mm Folded Size: 180\*170\*590mm




# Chair Stretcher Spine Board





SKB1C08 SKC-2B3 Chair Stretcher SKB1C09 SKC-2B3 Chair Stretcher

Model			
SKB1C08		1060*520*140mm	159kg
SKB1C09	1880*500*200mm	1370*500*200mm	159kg

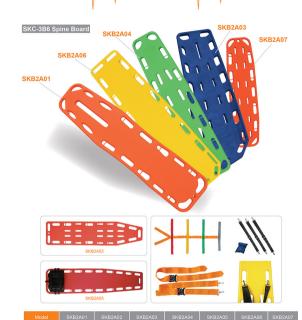



# Scoop Stretcher Basket Stretcher



#### SKC-3B7 Scoop Stretcher

#### Technical features:

- Equipped with a clutch hinge in the middle of both end, which can separate the stretcher into two parts.
- The patients can be quickly & easily shovelled on the stretcher and als stretcher can be taken out from the back of patients without moving the Can be used together with head immobilizer.
- Allow X-ray examination.
- · Suitable for patients with suspected hip, pelvic, and spinal injuries.
- Shell: High density polyethylene.

so the se patients.	C. T. T.


Model	SKB2B01	SKB2B02	SKB2B03	SKB2B04	SKB2B05
	2140*420*70	1850*440*60	1620*445*70	1720*430*70	1870*420*60
Folded Size(mm)	1670*420*70	,	1200*445*70	1230*430*70	/



SKC-3B7 Scoop Stretcher



**SKB2B03** SKC-3B7 Scoop Stretcher

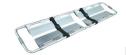

SKB2B02 SKC-3B7 Scoop Stretcher



**SKB2B04** SKC-3B7 Scoop Stretcher

#### SKC-3B10 Helicopter Stretcher/Basket Stretcher









### SKB2C02

	•	
Folded Size(L*W*H)	Load Bearing	N.W.
/	272kg	18kg

SKB2C01	2160*610*190mm	/	272kg	18kg
SKB2C02			272kg	





Size: 1890\*640\*420mm 25\*25mm square tube/25\*35 aluminum tube



SK-TB003 Foldable Tent Bed Size: 1950\*680\*425mm

20\*20mm square tube/1mm wall thickness



### Accessories



### SKB2D001 Cervical Collar

Size: L(630mm), M(570mm), S(510mm)



Cervical Collar



Cervical Collar



SKB2D002 Head Immobilizer

(A)Unfolded Size(L\*W\*H): 420\*280\*180mm (B)Unfolded Size(L\*W\*H): 340\*260\*140mm



SKB2D007 Cervical Collar

ength is adjustable





#### SKB2D102 Aluminum Mouldable Splint

Full Size(L\*W\*H): 1850\*470\*65mm Folded Size(L\*W\*H): 900\*1100mm



### Transfer Stretcher

Product Size(L\*W\*H): 1800\*560\*50mm Folding Size(L\*W\*H): 900\*560\*50mm





#### SKB2D201 Emergency Foil Blanket

Size: 1300\*2100mm 1400\*2100mm 1600\*2100mm



#### SKB2A08 CPR Board

Product Size(L\*W\*H); 640\*430\*80mm



#### SKB2A09 Rescue Tube

Product Size(L\*W\*H): 700\*240\*145mm



### Rescue Tube

Product Size(L\*W\*H): 1280\*160\*90mm 980\*160\*90mm



### Accessories Special stretchers



#### SKC-3D5 Traction Splint Set

- Relieve the tissue pressure on fractured leg avoiding further hemorrhage, shock and secondary nervous, vascular, bone and muscular damages. Made of aluminium alloy.
- Adjustable length, and nylon leg straps automatic traction
- Set of adult and pediatric sizes together.
- ◆ Product Size(L\*W\*H): 1380\*180\*170mm/1230\*180\*170mm
- Folded Size(L\*W\*H): 930\*70\*180mm/860\*70\*180mm Packing Size(1pc/carton): 950\*80\*220mm/880\*80\*220mm

### SKC-3C4 Body Splint

- Forehead strap provides support for fixing the head. The soft strap provides ergonomic support for fixing the chin.
- Provides both immobilization and comfort with a built-in plywood.
- spine board along the back. It is applied to rib injured.
- Three color-coded body straps with quick-clip, snap-lock buckles make it simple to apply in low light or cramped areas.
- . When used in conjunction with a cervical collar, it provides to a long spine board.
- Comes complete with two reusable, plastic-coated, head/chin straps. neck pad and durable carrying case.





### SKC-3C4 Pediatric Immobilization Strap

- Designed to adapt an ambulance cot to securely hold a child from 7 to 40lb(3.2 to18kg)
- Uses three restraint straps. · Fully-adjustable, five-point harness system fits children
- firmly and provides safe restraint during transport. · Rolls up compactly for efficient storage and retrieval.
- Made from nontoxic, easy-to-clean vinyl,
- Width(open): 480mm. Length(open): 800mm.
- · Rolled(diameter): 100mm, Rolled(length): 480mm.
- . Weight: 1kg.
- Load(min): 5kg. Load(max): 18kg.



Full Size(L\*W\*H): 2120\*750\*70mm Folded Size(L\*W\*H); 700\*600\*370mm







### SKC-3C4 Neil Robertson Stretcher

Hame: Hard wood Outer Material: Canvas Horizontal Safety: 120kg Vertical Safety: 120kg



- · Forehead strap provides support for fixing the head. The saft strap provides ergonomic support for fixing the chin.
- This product is applied to child's shoulder, chest, abdomen and leg's fixing, makes the child 's heart comfortable and respiration smooth.
- The fastened leg and ankle straps allow for separate treatment and immobilization of leg. Can be radiolucent for X-ray. Serves as a spine immobilizer during transport as well as a
- restraint system for children from 28" to 54" tall, weighting 20lb to 90lb. · Sewn-in lifting Handles at both ends for easy handling in confined areas.



### Isolation Stretcher Soft Stretcher



SKB3A008-1 is a complete shut-off negative pressure compartment, used for temporary isolation and transferring of respiratory tract infectious disease patient or suspected patient.

It can prevent pathogens from further spread; protect the medical staff and public environment. The fresh air is delivered to patient's cabin after through appropriate filter timely, provide comfortable security escort environment for patient.

#### Technical indicators and configuration:

The air filtration efficiency to 0.3 µm aerosol particles is not less than 99.99%. The effective time is 6 hours.

- Negative pressure device for opening and closing seal of the cabin, equipped with air filtration system and negative pressure generation system. The polluted air in cabin will be explicately fain after the filtration from canister and the negative pressure is established by adjusting the fair speed, so that to prevent pathogens from teating.
  Use lithium battery or arb tattery, how kinds of power-supply modes can be switched. Lithium battery can run continuously for
- 2 hours after fully charged.
  Suitable for use inside the ambulance, transferring infection patient.
- ◆ Cabin negative pressure difference inside and outside is not less than 40pa, negative pressure build time is less than 2 min.
- When negative pressure device at work, the noise of the cabin is less than 72 db (A)
   Ventilation volume more than 50 L/min

#### Accessories:

Cabin \*1 Canister \*1 Dust filtering box \*1(already installed in cabin) Air-way tube \*2(1pc long + 1pc short)
Fan \*1 Controller \*1 Lithium battery \*1 Lithium battery charger \*1 Car ignition extension cord \*1
Let a gloves \*1 Gloves ring \*1 User and maintenance manual \*1







Full Size(L\*W): 1870\*420mm Folded Size(L\*W): 880\*240mm



SKB3A102 SKC-3E5 Soft St

Full Size(L\*W): 1980\*970mm



Various colours of nylon ties for your choice, along with the safety belt.

	retcher ca packaged	
	, a a mag a a	aage

Model	Full Size(L*W)	Load Bearing	Heat Resistance	N.W.
SKB3A103	2250*500mm	160kg	-20°C~45°C	5.2kg

### First Aid Kit Resuscitators/Face Mask



SKB-5A002 First Aid Kit

Antibacterial bandage \*1 Medical tape \*1 Anriseptic cleansing wipe \*3





SKB-5B005 First Aid Kit







SKB-5A006 First Aid Kit

SK-EH021 First Aid Kit



SKB-5A004

First Aid Kit

SKB-5B007 First Aid Kit



SKB-5B020 First Aid Kit



SKB-5B010 First Aid Kit



SKB-5B012 First Aid Kit

# Silicone Manual Resuscitators

SKB-5C001 for adult and child SKB-5C002 for child and infant SKB-5C003 for infant and neonate

Reservoir bag Optional Package: PP box/Paper box/Poly bag



### Disposable PVC Manual Resuscitators

SKB-5C006 for adult and child SKB-5C007 for child and infant SKB-5C008 for infant and neonate

Accessories: PVC face mask, Oxygen tubing, Reservoir bag Optional Package: PP box/Paper box/Poly bag Optional Color: Green/ Blue







### SKR-5C010 without check valve



#### SKR-5C015 PVC Face Mask(0#-5#) Guedel Airway(4#-12#)





#### SEBS Manual Resuscitators

SKB-5C011 for adult and child SKB-5C012 for child and infant SKB-5C013 for infant and neonate

Accessories: PVC face mask. Oxygen tubing. Reservoir bag Optional Package: PP box/Paper box/Poly bag Optional Color: Green/ Blue





#### SKB-5C014 **CPR Face Mask**

CPR face mask/Packaged by poly bag CPR face mask/Packaged by PP case Optional Accessories: Alcohol pad. PVC glove Optional PP Case Colour: White/Red/Green/Orange/Blue







# Oxygen Intake Device Tent Bed



SK-EH041 Oxygen Bag





Oxygen Cylinder

Material: 37Mn Steel Weight: 2.5kg



SK-EH005 Oxygen Cylinder



SK-D001 Defibrillator



SK-D003 Defibrillator







Medical Oxygen Flowmeter Medical Oxygen Flowmeter Medical Oxygen Flowmeter Medical Oxygen Flowmeter



SK-EH035





Field Camp

For: 3-4 People Style: 2 People, 3 People, 5 People Size: Customized Pole Material: Aluminum



Relief Tent

Type: First aid tent Size: Customized



Flow Meter



SK-EH048 Flow Meter



SK-EH049



SK-EH050 Medical Oxygen Flowmeter Argon Gas Pressure Reducer



SKB-4A003 Tent



SKB-4A007 Tent




SKB-4A008 Tent



SKB-4B005 Tent

# Emergency Traning Model Mortuary Refrigerator

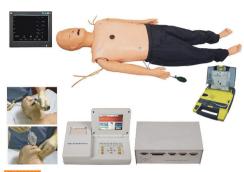

### SKB-6A002

Comprehensive Emergency Skills Training Manikin

- Full body manikin \*1
- ACLS monitor \*1 ACLS converter \*1
- ◆Aluminum & plastic box \*1
- . CPR operation pad \*1 ◆ Exchangeable lung bag \*5
- CPR face shield sheet \*50pcs/box
- Temperature sensing print paper \*2 • KAR/S7 BP training arm \*1





SKB-6A012 Child Trachea Intubation Model



SKB-6A003 CPR Training Manikin



SKB-6A004



SKB-6A007



SKB-6A008

Half Body CPR Training Manikin Infant CPR Manikin Adult Obstruction Model



Mortuary Refrigerator



SKB-7A003 Mortuary Refrigerator

	Wortdary Reingerator			mortuary itemigerator			
Model		SKB-7A002		SKB- 7A004	SKB-7A005		
Lowest Temperature	-10~-30°C	1~5℃	1~5℃	1~5℃	1~5°C	-10~5℃	
Quantity of Compressorunit	30°C	-16~-30°C	-16~-30°C	-16~-30°C	-16~-30°C	1	
Compressor Ploidy	1 set	2 sets	3 sets	4 sets	6 sets	1	
Computer					6 sets		
Power Supply	Individually, AC220V 50HZ						
Compressor Type	Individually, AC220V 50HZ / Imported closed /						
Corpse Tray Size(mm)	1900*600*360(Depth, Width, Height)						
Dimension(mm)		2500*850*1210		2500*1600*1210	2500*1600*1640		
Total Weight	110kg	230kg	320kg	500kg	600kg	102kg	
Relative Humidity		1		,	/		







SKB-7A002



SKB-7A005



SKB-7A006 Mortuary Refrigerator

# Body Bag Church Trolley







SKB-7B002 Body Bag



SKB-7B003 Body Bag



SKB-7B002(2) Body Bag

		Locking Positio	E. 11 6			
Model	1st position	2nd position	3rd position	Folding Size(L*W*H)	Bearing Load	
	710*550*720mm	1260*550*630mm	1550*550*560mm	560*210*730mm	450ka	









SKB-7C003

Church Trolley(Aluminum Alloy)



#### SKB-7C006 Urn Holder

Detailed product description: It is made of steel material. Spray-coated.

Spray-coated.

It is mainly used for transferring te cinerary casket.

Commonly used to fix the cinerary casket in
the car during transportation.



#### SKB-7C008 Coffin Trolley

The coffin rack is made of high strength aluminum alloy material, it is mainly used for coffin placing.



#### SKB-7C007 Church Trolley

Technical characteristics: Made of solid durable steel, foldable and portable framework,

easy to operate, the length can be adjusted freely, two small wheels with brake can be used alone.

#### Specifications: Unfold Size(L\*W\*H): 1220\*50\*55\* Packing Size(L\*W\*H): 1280\*520\*

Packing Size(L\*W\*H): 1260\*520\*70mi N.W.: 13.5kg G.W.: 15.5kg Load Bearing: s450kg



#### SKB-7C009 Catafalque

It is made of steel pipes, three level plating-coated, it is also easy for cleaning. It can save a lot of space when it is folded, and it is mainly used for coffin placing.

SKB-7C008	1150*550*630mm	1	≤159kg	10kg	13kg
			≤350kg		