



CR0q

Children Bed

SK-A2



SAIKANG / Medical

LEADING BRAND OF CHINA MEDICAL FURNITURE

Stock Code: 870098

TECHNICAL PARAMETERS

External L*W	1600*740*1040mm
Bed Height	500mm
Mattress Platform	1580*660mm

TECHNICAL CONFIGURATION

5" Covered Castors	4pcs
Full Length Side Rail	1 set
Integrated Bed Ends	1 set
IV Pole	1 pc
IV Pole Prevision	4 pcs
Drainage Hook	2 pcs

STANDARD ACCESSORY

One telescopic IV pole

OPTIONAL ACCESSORY

Basic foam mattress - 12cm thickness
Anti-decubitus foam mattress - 14cm thickness
Vertical oxygen holder
Single arm lifting pole
Utility Shelf

QUALITY CONTROL



RICH PARTS

LEADING BRAND OF
CHINA MEDICAL FURNITURE



Punching platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform



Powder coated steel frame, MDF material fire-proof plate. size of the head&foot board can be customized. With carton photograph




Foldable guardrail with protection handle, modern and beautiful, easy to install, easy to clean, anti-pinch, protect patient's safety



TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time



Full SS IV pole, bearing 15kg

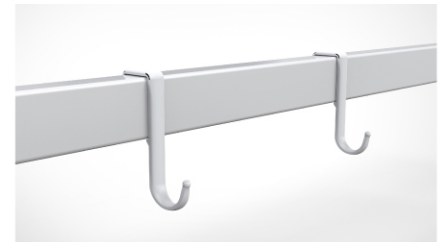

Plastic drainage hook, movable
Shoes holder: 800*400mm, adjustable
Oxygen cylinder rack - adjustable for anysize tank
Lifting pole - Bearing 765N force





INNOVATION OF SAIKANG

LEADING BRAND OF CHINA MEDICAL FURNITURE



Design Innovation



Research Innovation



Manufacturing Innovation



Maintenance Innovation



Service Innovation

JIANGSU SAIKANG MEDICAL EQUIPMENT CO.,LTD.

24hrs: 0086-18021231901 
www.chinasaikang.com

SAIKANG (China Headquarter)

No.35 Lehong Road,Zhangjiagang City,
Jiangsu Province,China.
TEL: +86-512-58525850
FAX: +86-512-58520713
E-mail: inquiry@jssaikang.com

SAIKANG (USA Branch)

Saikang Medical Technology, Inc
17800 Castleton St Ste 180,
City of Industry, Ca 91748
Tel: +1 949648-2776

SAIKANG (UK Branch)

Uk Healthward International Limited
Owner(s) Address DEPT 114,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +44 2032399738

SAIKANG (East Africa Branch)

Saikang Medical Solutions (East Africa) Ltd.
Viraj Complex Unit 3 Mombasa Road,
PO Box 14262-00800, Nairobi.
TEL: +254 733-704 690

SAIKANG (Mid East Branch)

Saikang Medical Supplies (Mid East) LLC
X22-S03-England Cluster,
International City, Dubai UAE
TEL: +20 128 8962497



Contact



Website

As the circumstance of printing, the photos might have minor difference on color, tactile from real products.
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.