



H6k

ELECTRIC HOMECARE BED

Registration code: SKD-C




SAIKANG / Medical

Healthier
your life

H6k

FEATURES AND HIGHLIGHTS

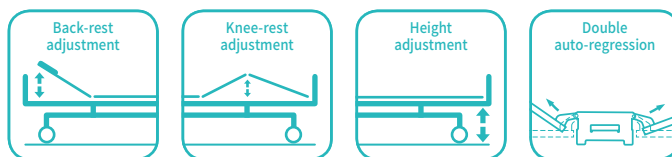

Technical parameters

▪ External size (LxWxH)	2015x1050x360-800 mm
▪ Mattress platform	1925x820 mm
▪ Safe working load	200 kg
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)
▪ Voltage	100-240 V
▪ Frequency	50/60 Hz
▪ Plug as required (hospital grade)	

Technical configuration

<input checked="" type="checkbox"/> Linear Motor	3pcs
<input checked="" type="checkbox"/> 5" Covered Castors	4pcs
<input checked="" type="checkbox"/> Foldable Side Rail	1set
<input checked="" type="checkbox"/> Wood Bed Ends	1set
<input checked="" type="checkbox"/> Wood Bed Side Panel	1set
<input checked="" type="checkbox"/> IV Pole	1pc
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs
<input checked="" type="checkbox"/> Mattress Holder	2pcs



Steels | Painting

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50kg impact.

**Welding | Process**

It adopts imported intelligent robot welding and forming, 360° welding without dead angle.

**Electric System**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Auto-Regression**

When repositioning the backrest and thigh rest, the patented design of the Auto-Regression of the bed surface minimizes the pressure and shear force on the pelvic area of the patient, this results in a significant reduction in the development of pressure ulcers.

**Wooden Head & Foot Board**

The headboard and footboard of solid wood is used for firm fastening, good load-bearing ability.

**Platform**

Stamping forming platform, smooth surface no welding scar, bearing 200KG load. Option resin X-ray translucent platform.

**Guardrail**

It adopts aluminium alloy foldable guardrail, The guardrail switch is made of new high-strength ABS engineering plastic injection molding.

**Castors**

5" covered casters. TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts. Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.

**Drainage Hooks**

Laser cutting drainage hook, both sides.

**IV Pole**

Full SS IV pole, bearing 15kg.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue