




# SKB038-4

Patient  
Trolley



**SAIKANG** / Medical

LEADING BRAND OF CHINA MEDICAL FURNITURE  
Stock Code: 870098

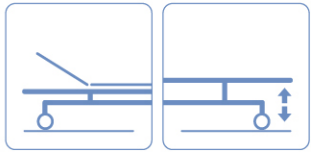

Back-rest

Height Adjustment

## MAIN TECHNICAL PARAMETERS

Size	1960*700*(620-880)mm
Back-rest	0-75°
Height Adjustment	620-880mm

## TECHNICAL CONFIGURATION

Crank System	1 set
Back rest Air spring	1 pc
5'covered casters	4 pcs
Sixth-foldable aluminium alloy guardrail	2 pcs
bumper straps	1 set
6cm PU leather mattress	1 pc
stainless steel I.V pole with 2 hooks	1 pc
Oxygen cylinder rack (Optional)	1 pc

## QUALITY CONTROL

### STEELS | PAINTING

Baosteel from Fortune Global 500  
I I Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact

0.12mm | 60° | 50kg

### WELDING | CHASSIS

Panasonic Robotic ensure, 360° full smooth welding  
Solid bed base made of 5.0mm stamped rotating plate  
Patent angle indicator

360° | 5.0mm

### BUMPER STRAPS

The trolley frame with the bumper straps could reduce the impact force

### TROLLEY SURFACE

The trolley surface is made of powder coated steel, smooth surface no welding scar

## RICH PARTS


LEADING BRAND OF  
CHINA MEDICAL FURNITURE




Steel hanging crank, 100,000 times remaining intact, it can be fold under the bed when not in use



5 inch  
TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time

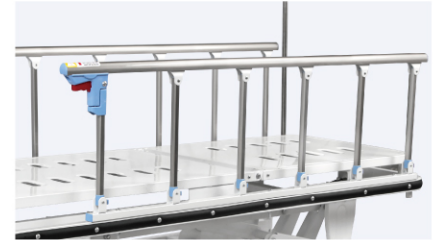

Aluminum alloy side rail, 6-staff, collapsed type with safe lock

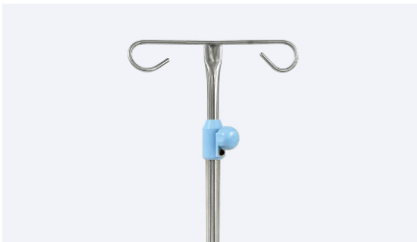

Adjustable stainless steel I.V. pole, outside diameter  $\Phi 19*1$ , inside diameter  $\Phi 16*1$ , telescoping scope 750mm, load capacity: 15KG



6cm thickness flame retardant high-density foam with U leather cover



Oxygen cylinder rack - adjustable for any size tank





## INNOVATION OF SAIKANG

LEADING BRAND OF CHINA MEDICAL FURNITURE



Design Innovation



Research Innovation



Manufacturing Innovation



Maintenance Innovation



Service Innovation

### JIANGSU SAIKANG MEDICAL EQUIPMENT CO.,LTD.

24hrs: 0086-18021231901   
[www.chinasaikang.com](http://www.chinasaikang.com)

#### SAIKANG (China Headquarter)

No.35 Lehong Road,Zhangjiagang City,  
Jiangsu Province,China.  
TEL: +86-512-58525850  
FAX: +86-512-58520713  
E-mail: [inquiry@jssaikang.com](mailto:inquiry@jssaikang.com)

#### SAIKANG (USA Branch)

Saikang Medical Technology, Inc  
17800 Castleton St Ste 180,  
City of Industry, Ca 91748  
Tel: +1 949648-2776

#### SAIKANG (UK Branch)

UK Healthward International Limited  
Owner(s) Address-DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +44 2032399738

#### SAIKANG (East Africa Branch)

Saikang Medical Solutions (East Africa) Ltd.  
Viraj Complex Unit 3 Mombasa Road,  
PO Box 14262-00800, Nairobi.  
TEL: +254 733-704 690

#### SAIKANG (Mid East Branch)

Saikang Medical Supplies (Mid East) LLC  
X22-S03 England Cluster,  
International City, Dubai UAE  
TEL: +20 128 8962497



Catalogue



Website

As the circumstance of printing, the photos might have minor difference on color, tactile from real products.  
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.