



A202

ELECTRIC
OPERATION TABLE

Registration code: SKL-A




SAIKANG / Medical

Healthier
your life

A202

FEATURES AND HIGHLIGHTS

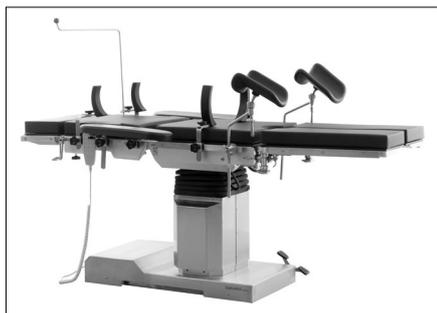

Technical parameters


External size (LxWxH):	2000x500x(750-1050)mm
Reverse trendelenburg:	0-22°
Trendelenburg:	0-25°
Lateral tilt left/right:	0-20°
Head plate up/down:	45°/90°
Back plate up/down:	70°/35°
Leg plate down:	0-90°
Leg plate outwards:	0-90°
Kidney bridge elevation:	0-80mm
Horizontal sliding:	0-300mm
Safe working load:	360kg

Technical configuration

<input checked="" type="checkbox"/> Arm support	2 pcs
<input checked="" type="checkbox"/> Leg support	2 pcs
<input checked="" type="checkbox"/> Body Support	2 pcs
<input checked="" type="checkbox"/> Shoulder Support	2 pcs
<input checked="" type="checkbox"/> Slider	7 pcs
<input checked="" type="checkbox"/> Tightener	2 pcs
<input checked="" type="checkbox"/> Anaesthetic Screen	1 pc
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Handwheel	1pc

Material

The chassis cover is made of high quality Ni-Cr stainless alloy steel by seamless robot welding, corrosion resistance, disinfection resistance, rust-proof and durable.

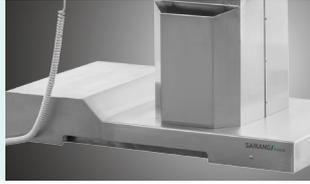

Mattress

Removable memory mattress, good rebound resilience, anti-static electricity and anti-microbial.



Control System

Electric-hydraulic control system, non-gear structure, all functions are controlled by the hand controller.



Braking system

Mechanical brake device at the side of the base, with fixed or moving function, flexible and reliable.



Hand controller

All functions are controlled by the hand controller.



Kidney bridge

Kidney bridge is adjusted by the hand wheel at the side of head board, more convenient, with separate mattress.



Battery could support

Battery could support about 50-80 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.



Horizontal movement

300mm horizontal movement function, with high strength phenolic resin table top, could be used for X-ray and C-arm.



Leg plate

The up and down folding of the leg plate is manually controlled by the high strength gear, which makes the operation more effort-saving and smooth.



Armrest



Anesthesia Frame



Leg Crutches



Shoulder Support



Body Supporter



Professional functions

Lowest position of table top

750mm



Highest position of table top

1050mm



Back plate up

0°-70°



Back plate down

0°-35°



Reversed trendelenburg

0°-22°



Trendelenburg

0°-25°



Right lateral tilt

0°-20°



Left lateral tilt

0°-20°



Head plate upward

0°-45°



Head plate downward

0°-90°



Leg section evagination

0°-90°



Leg plate up/down

-45°-90°



Kidney bridge elevation

0-80mm



Horizontal movement

0-300mm





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue