



A99-13

ELECTRIC OBSTETRIC BED

Registration code: SKF-C




SAIKANG / Medical

Healthier
your life

A99-13

FEATURES AND HIGHLIGHTS

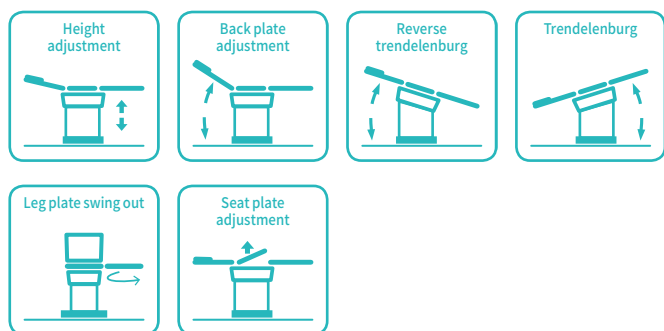

Technical parameters

▪ External (LxW):	2100x600mm
▪ Table height:	750-1000mm
▪ Trendelenburg:	0-22°
▪ Reverse trendelenburg:	0-22°
▪ Back plate up:	0-60°
▪ Seat plate up:	0-16°
▪ Leg board outwards:	0-90°
▪ Arm board outwards:	0-180°
▪ Load bearing:	≤300kg



Technical configuration

<input checked="" type="checkbox"/> Armrest	2 pcs
<input checked="" type="checkbox"/> Knee Curtch	2 pcs
<input checked="" type="checkbox"/> Pull handle	2 pcs
<input checked="" type="checkbox"/> Anesthesia frame	1 pc
<input checked="" type="checkbox"/> SST rinsing basin	1 pc
<input checked="" type="checkbox"/> Pedal controller	1 pc
<input checked="" type="checkbox"/> LED examination light	1 pc



Material

Base structure is made of high-quality 304 stainless steel, robot welded with high strength and modern appearance.

**Mattress**

Polyurethane foam mattress, waterproof and seamless, easy cleaning and disinfection.

**Foot pedal controller**

All movements are controlled by foot pedal controller.

**Brake**


Pedal brake, easy to move and lock. Pedal down to release the inner castors down. Pedal up to retract the castors, which enables the table locked on the floor.

**LED examination light**

Standard accessory: LED examination light. More convenient to do gynecological examination.

**Arm plate**

Flexible arm plate adjustment.

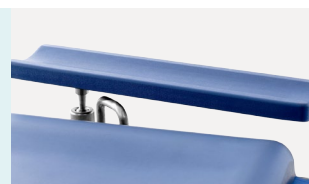**Battery could support**

Battery could support about 80-120 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table runing in the event of power failure.

**Leg supports**

Equipped with two different leg supports to meet different surgical needs.

**Auxiliary plate**

Auxiliary plate is easy for rotating and dismantling. Rotating angle: 90°.

**Pull handles**

Pull handles can be folded to horizontal position. Easy patient mobilisation.

**Sewage basin**

Equipped with stainless steel sewage basin and other various accessories.

**Electrical controlled**

All movements are electrical controlled: Back plate, seat plate, table elevation, trendelenburg and reverse-trendelenburg.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue