



# B8e

## MULTIFUNCTIONAL ICU ELECTRIC BED

Registration code: SKD-C






**SAIKANG**/Medical

Healthier  
your life

## B8e

### FEATURES AND HIGHLIGHTS

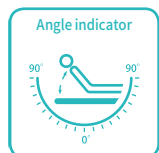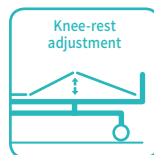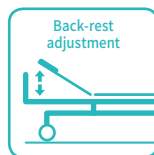

### Technical parameters

|                                     |                      |
|-------------------------------------|----------------------|
| ▪ External size (LxW):              | 2170x1065 mm         |
| ▪ Mattress platform:                | 1925x900 mm          |
| ▪ Safe working load:                | 250 kg               |
| ▪ Height adjustment:                | 435-725 mm           |
| ▪ Back-rest adjustment:             | 0-75° (±10°)         |
| ▪ Knee-rest adjustment:             | 0-35° (±10°)         |
| ▪ Reverse Trendelenburg:            | 0-12°                |
| ▪ Trendelenburg:                    | 0-12°                |
| ▪ Angle indicator:                  | 0-90°                |
| ▪ Voltage and frequency:            | 100-240 V · 50/60 Hz |
| ▪ Plug as required (hospital grade) |                      |




### Technical configuration

|                                                                                |       |
|--------------------------------------------------------------------------------|-------|
| <input checked="" type="checkbox"/> Linear Motor                               | 4 pcs |
| <input checked="" type="checkbox"/> 5" Double Side Castors, Central Locking    | 4 pcs |
| <input checked="" type="checkbox"/> ABS Side Rail Extra Height                 | 1 set |
| <input checked="" type="checkbox"/> ABS Bed Ends with Integrated Control Panel | 1 set |
| <input checked="" type="checkbox"/> IV Pole                                    | 1 set |
| <input checked="" type="checkbox"/> IV Pole Prevision                          | 4 pcs |
| <input checked="" type="checkbox"/> Drainage Hook                              | 2 pcs |
| <input checked="" type="checkbox"/> Central Locking Pedal                      | 2 pcs |
| <input checked="" type="checkbox"/> Electric CPR                               |       |




### Optional

|                                                |  |
|------------------------------------------------|--|
| <input type="checkbox"/> Manual CPR            |  |
| <input type="checkbox"/> Weighing System       |  |
| <input type="checkbox"/> Back-up Battery       |  |
| <input type="checkbox"/> X-ray Bed Surface     |  |
| <input type="checkbox"/> Monkey Pole           |  |
| <input type="checkbox"/> Balkan Traction Frame |  |
| <input type="checkbox"/> Directional Wheel     |  |

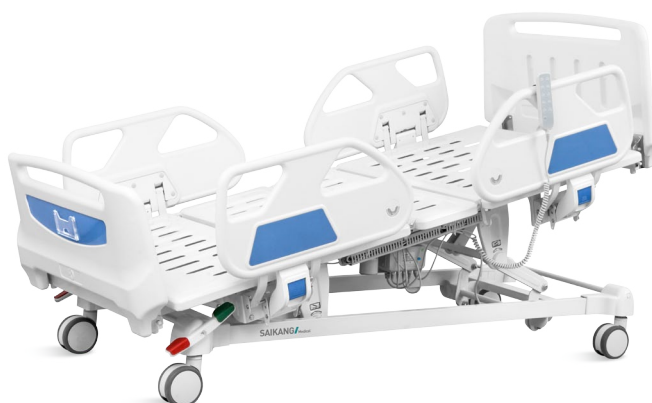

**Design**

E-shaped large guardrail, with angled display and built-in controller, meets safety standards, providing comforting feeling for patients.

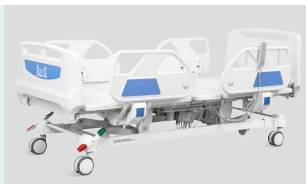**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding. Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**


Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Virgin ABS**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**Angle Indicator**

Patented angle indicator.

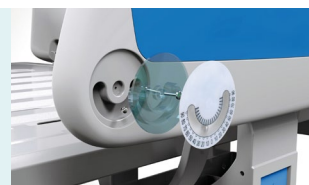**Wheels**

Secure® Castors  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**Durable Brakes**

Break spring keep qualified after 10,000 times uninterrupted break.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Side Rail**

Full ABS material, tuck away side rail with built-in controller. From top of the guardrail to bed surface is about 345mm, more height provides better protection for patient. Bottom can endure 100,000 times press, buffer, built-in angle indicator.

**Head & Foot Board**

Full ABS material, Bottom can endure 100,000 times press  
Head and foot board can be reversible.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic hook.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SAIKANG** Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎ +86 18021231901 · 🕒 24/7  
✉ sales@chinasaikang.com  
🌐 www.saikangmedical.com



Company video



General catalogue