



CD3q

THREE FUNCTIONS CHILDREN BED

Registration code: SK-CD1




SAIKANG / Medical

Healthier
your life

CD3q

FEATURES AND HIGHLIGHTS

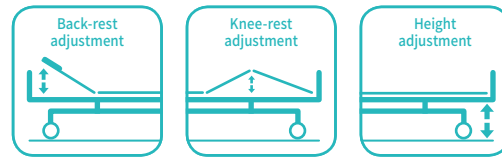

Technical parameters

▪ External size (LxW)	2030x1045 mm
▪ Bed height	500-770 mm
▪ Mattress platform	1705x900 mm
▪ Back-rest adjustment	0-75° ($\pm 10^\circ$)
▪ Knee-rest adjustment	0-35° ($\pm 10^\circ$)



Technical configuration

<input checked="" type="checkbox"/> Manual Cranks	3pcs
<input checked="" type="checkbox"/> 5" Hardshell Castors	4pcs
<input checked="" type="checkbox"/> Aluminium Alloy Full Length Side Rail	1set
<input checked="" type="checkbox"/> ABS Grace Bed Ends with Safe Lock	1set
<input checked="" type="checkbox"/> IV Pole	1pc
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs



Design

Head foot board, with curve surrounded design, provide comforting feeling for patients.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding Solid bed base made of 5.0 mm stamped rotating plate.

**Hand Control**

Retractable hand crank, 100,000 times remaining intact Extreme protection device, sounds when tipped to the max position.

**Platform**

Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.

**Guardrail**

Sliding full side rails with distinctive design, elegant appearance; The lock bar is in conspicuous red, easy operation, sliding upper and down smoothly; Anti-pinch design, effectively protect the safety of children in use.

**Head & Foot Board**

Head and foot board can be reversible, extra height design. ABS Grace Bed Ends with Safe Lock, automatic lock without any manual operation on lock, on/off removable mark on hook lock.

**Castors**

5" hardshell casters
TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts.
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic hook.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue