



SKS025

PHENOLIC RESIN
BEDSIDE TABLE



SKS025

FEATURES AND HIGHLIGHTS

Technical parameters

- External size (LxWxH) 475x485x810 mm

Technical configuration

<input checked="" type="checkbox"/> Top Guardrail	Stainless Steel
<input checked="" type="checkbox"/> Drawer	1pc
<input checked="" type="checkbox"/> Cabinet	1pc
<input checked="" type="checkbox"/> Casters	4pcs

Handle

The handle is made of aluminum alloy with surface anodized, which is firm, rust-proof and easy to clean.



Casters

Casters is using medical silent material, double side casters, stable and noiseless. 2 casters with brake, lockable at 360 degree, 2 casters free and swivel to any direction.




Frame

Aluminum alloy 4 column frame, stable and firm.



Material

Table body 3 side walls are made of powder coated steel, table top, drawer and door are made of phenolic compact board with rounding treatment, anti-acid, anti-alkali, waterproof and anti-seepage, easy to clean and flame retardant for 4 hours.

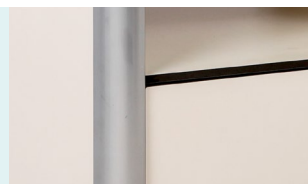

Structure

A drawer, a cabinet, a storage open case in the middle.



Sliding Track

The sliding track adopts silent damping track, which can effectively prevent the drawer from sliding out in the course of moving the cart.



Color

Option: color (color samples) .





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue