



D6w

THREE FUNCTIONS ELECTRIC BED

Registration code: SKD-C




SAIKANG / Medical

Healthier
your life

D6w

FEATURES AND HIGHLIGHTS

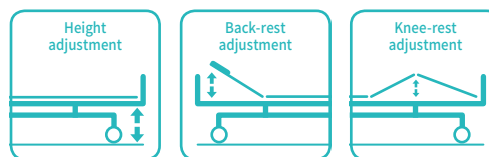

Technical parameters

▪ Size (LxW)	2245x1000 mm
▪ Bed height adjustment	435-745 mm
▪ Mattress platform	1925x830 mm
▪ Safe working load	250 kg
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)



Technical configuration

<input checked="" type="checkbox"/> Linear motor	3 pcs
<input checked="" type="checkbox"/> Hand controller	1 set
<input checked="" type="checkbox"/> 5" hardshell castors with break	4 pcs
<input checked="" type="checkbox"/> Foldable side rail with ProErgo-Handle™	1 set
<input checked="" type="checkbox"/> ABS bed ends with safe lock	1 set
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> IV pole prevision	4 pcs
<input checked="" type="checkbox"/> Drainage hook	2 pcs



Optional

<input type="checkbox"/> Manual CPR	
<input type="checkbox"/> Back-up battery	
<input type="checkbox"/> X-ray translucent mattress platform	
<input type="checkbox"/> Directional wheels function	
<input type="checkbox"/> Monkey pole	
<input type="checkbox"/> Aluminium balkan frame	



Design

Extra Long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**


Panasonic robotic ensure, 360° full smooth wedding. Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

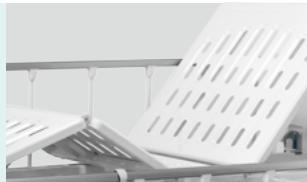**Side Rail**

L-shaped foldable guardrail (Patent), bended tube, solid installing, 0 shaking, anti-pinch, protect ensure patient's safety.

**Head & Foot Board**

Full ABS material, Bottom can endure 100,000 times press. Head and foot board can be reversible.

**Wheels**

TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts. Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**IV pole**

Stainless Steel Column with plastic hook.

**Drainage Hooks**

Mobile plastic drainage hooks.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue