



Oasis Series

# Y6t8y

## ELECTRIC BED

Registration code: SKD-C




**SAIKANG**/Medical

Healthier  
your life

# Y6t8y

## FEATURES AND HIGHLIGHTS

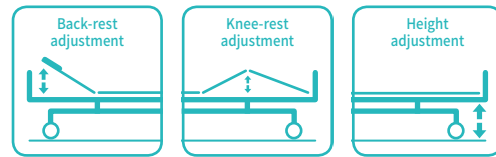

### Technical parameters

▪ External size (LxW)	2130x1060 mm
▪ Mattress platform	1925x900 mm
▪ Safe working load	250 kg
▪ Height adjustment	475-760 mm
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)
▪ Angle indicator	0-90°
▪ Position line	30°
▪ Voltage and frequency	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	
▪ Hand controller	

### Technical configuration

<input checked="" type="checkbox"/> Linear motor	3 pcs
<input checked="" type="checkbox"/> 5" double side castors, central locking system	4 pcs
<input checked="" type="checkbox"/> ABS Oasis side rail extra height	1 set
<input checked="" type="checkbox"/> ABS Oasis bed ends with safe lock	1 set
<input checked="" type="checkbox"/> Stamping forming bed platform	1 pc
<input checked="" type="checkbox"/> Plastic socket inside pp guardrail for phone or glasses	2 pcs
<input checked="" type="checkbox"/> Patent angle indicator for back-rest and trendelenburg position	4 pcs
<input checked="" type="checkbox"/> Ø70mm wall bumper wheel	4 pcs
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> IV pole prevision	4 pcs
<input checked="" type="checkbox"/> Drainage hook	2 pcs

### Optional

<input type="checkbox"/> Integrated control panel in side rail
<input type="checkbox"/> Back up battery
<input type="checkbox"/> Directional wheels function
<input type="checkbox"/> ABS mattress stopper
<input type="checkbox"/> Manual CPR
<input type="checkbox"/> Electric CPR



**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth wedding.  
Solid bed base made of 5.0 mm stamped rotating plate.

**Steels | Painting**

Baosteel from Fortune Global 500.  
11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Electric motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load.

**Wheels**

Castors with TPR tire, no worn-out after running 30KM, united forming without bolts.  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**Brake**

Three position braking system, transporting can maintain a straight line.

**Socket**

On side rail for placing the patient cellphone or glasses.

**Angle indicator**

Patent angle indicator, clearly show back-rest and tredelenburg position.

**Head & foot board**

Detachable PP head&foot board(Y19), with medical record card on foot panel, patent design, using the eccentric locking structure, one step quick removable.

**Siderail**


4 pcs PP siderail (T20), with patent angle indicators on each one, clearly show back-rest and tredelenburg position.

**Controller**

Hand controller within reach.  
Clear and easy to use controllers which is at the patient's disposal.

**Wall bumper wheel**

Φ70mm bumper wheels at the outermost periphery of the bed to protect the bed and wall, anti-collision and durable.

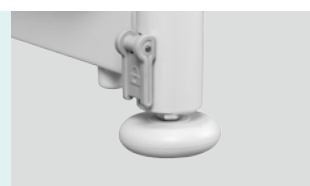**IV pole**

Stainless Steel Column with plastic hook.

**Drainage hooks**

Mobile plastic drainage hooks.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎ +86 18021231901 · 🕒 24/7  
✉ sales@chinasaikang.com  
🌐 www.saikangmedical.com



Company video



General catalogue