





# SKR-R05

## TABLET CART

Registration code: SKM-A

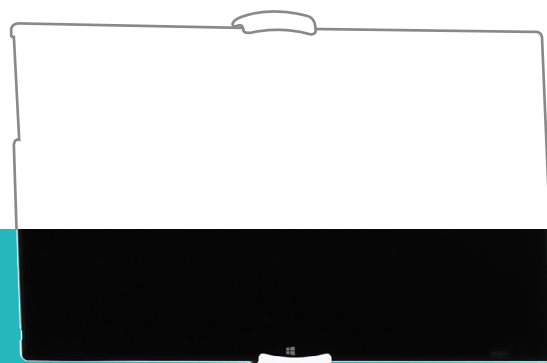

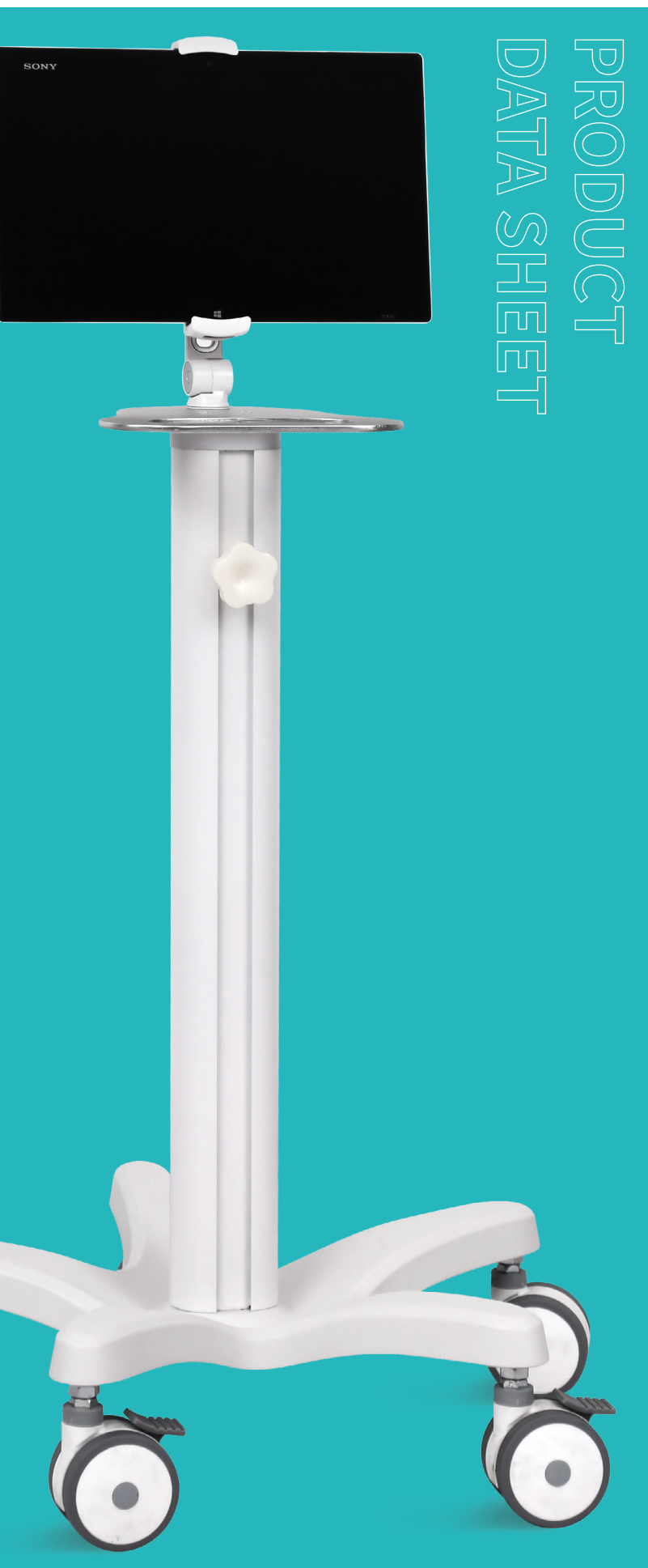

# SKR-R05

## FEATURES AND HIGHLIGHTS

### Product Introduction

The mobile doctor workstation based on wireless network integrates highly integrated computer, wireless communication, intelligent power management, mobile management technology, and is suitable for a variety of clinical work environments; whether it is outpatient, emergency, inpatient ward, operating room or ICU, it can help clinical Medical workers realize intelligent mobile ward rounds.

### Technical Parameters

▪ Dimensions W×H	510×830mm
▪ Work table L×W	210×200mm
▪ Overall weight	8.3kg
▪ Lifting method	manual lifting
▪ Lifting size	200mm
▪ Caster design	5 inch medical grade mute wheel, pedal brake



## Material

The table top and column are made of aviation-grade aluminum alloy, the base is made of one-time molding all-cast aluminum material, the whole car adopts high-voltage electrostatic spraying process, and imported environmentally friendly antibacterial plastic powder is used to ensure that the product meets the medical and health requirements, and the surface is matte and non-dazzling.




## Jaws

The tablet fixing clip can be rotated 360° and adopts a spring-type fixing method, which can be applied to tablet computers of different sizes. The buffer rubber inside the fixing clip can fix the tablet more stably and prevent scratches on the paint surface.



## Column

The column is equipped with a manual lifting function, which is easy to rise and fall, and can be quickly adjusted to a suitable height. It is locked by a knob, and can adapt to different sitting and standing postures, with a stroke of up to 20cm.

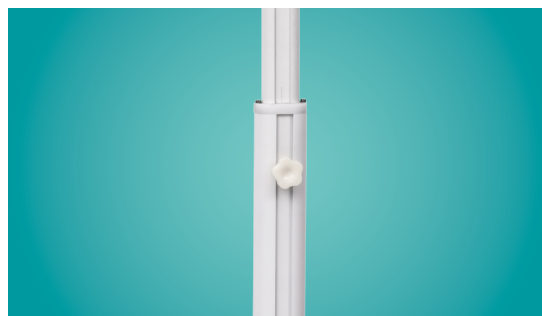

## Chassis | Casters

All-metal chassis, one-time cast aluminum molding, to ensure the center of gravity and level of the vehicle.

5-inch medical grade silent wheels, all with brakes, after tens of thousands of kilometers of long-term durability tests, flexible and durable.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### South America Branch

**Saikang Medical Technology, Inc**  
W3400 Corrientes, Argentina  
TEL: +55 (11) 2758-0438  
Rua Universal, Hollywood,  
São Bernardo do Campo,  
Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

**Saikang Medical Technology Ltd**  
Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

**Saikang Medical Technology Ltd**  
Sala Thammasop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

**Uk Healthward International Limited**  
DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 202111106

### Russia Branch

**LLC 'SAIKANG'**  
454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya Naberezhnaya,  
Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

**Saikang Medical Supplies (Mid East) LLC**  
Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

**SAIKANG** / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue