



X05-3

CHILDREN BED

Registration code: SK-C2




SAIKANG / Medical

Healthier
your life

X05-3

FEATURES AND HIGHLIGHTS

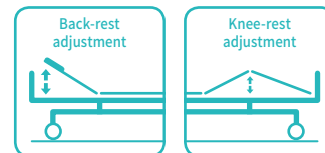

Technical parameters

▪ External size (LxW):	1800x800 mm
▪ Bed height:	700 mm
▪ External height:	1100 mm
▪ Mattress platform:	1435x690 mm
▪ Back-rest adjustment:	0-70° (±10°)
▪ Knee-rest adjustment:	0-40° (±10°)

Technical configuration

<input checked="" type="checkbox"/> Manual crank	2 pcs
<input checked="" type="checkbox"/> 5" covered castors	4 pcs
<input checked="" type="checkbox"/> Stainless steel full length side rail	1 set
<input checked="" type="checkbox"/> Stainless steel bed ends	1 set
<input checked="" type="checkbox"/> IV pole	1 pc
<input checked="" type="checkbox"/> IV pole prevision	4 pcs
<input checked="" type="checkbox"/> Drainage hook	2 pcs


Material

Stainless steel frame, baosteel from Fortune Global 500.




Bed platform

Mattress Platform: 1435*690mm
Stainless steel material, easy to clean and sterilize.

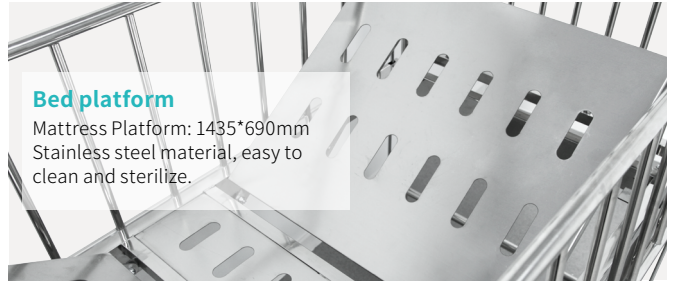


Guardrail

The stainless steel full-length guardrail protects the safety of children and patients in an all-round way. The guardrail can be raised and lowered, making it easy for patients to get out of bed and nurses to care.

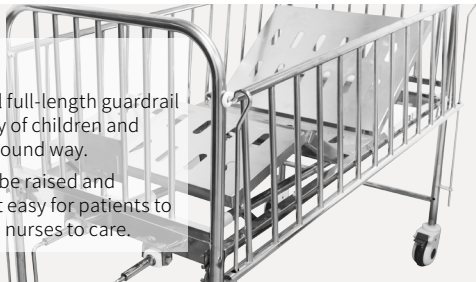

Crank

The retractable handle can be hidden to save space and prevent bumps.

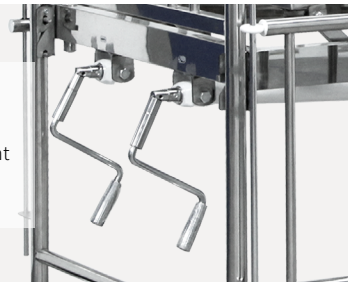

Castors

5" covered castors, two with brakes, easy to install, stable and smooth.



IV pole

Infusion rod with four hooks, bearing 10kg





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd
Sala Thammasop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited
DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya Naberezhnaya,
Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue