



SKE001-5

ACCOMPANY CHAIR



SAIKANG / Medical

Healthier
your life

SKE001-5

FEATURES AND HIGHLIGHTS

Technical parameters

- External size 650×2060×450mm
- Seat size 780×650×450mm(backrest height: 1060mm)
- Safe working load 200kg
- Color MOQ>20 colors can be customized



Design

Ergonomic design, stylish and comfortable; the chair can be adjusted as needed, convenient to attend the needs of patient, can be used as a sitting chair.




Frame

Reinforced design, high bearing capacity, the outer frame is made of bent and welded 38*1.2 round tubes, the chair legs are made of bent and welded 25*1.2 round tubes, safe load capacity up to 200kg.



Surface material

PU leather, moisture and water proof, thickened sponge inside, comfortable and high wear resistance.

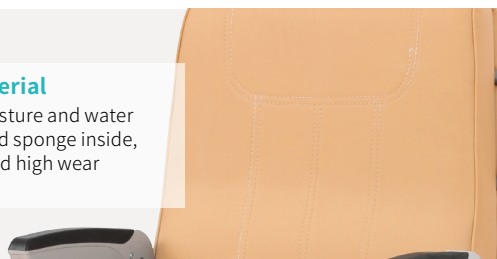

Armrest

Adjustable and foldable aluminum alloy armrest, PU leather surface.




Caster

Equipped with 3-inch silent caster*2pcs, 1.5-inch silent caster*4pcs, smooth, easy operable, the chair feet without casters are equipped with plastic anti-slip feet cover to protect the floor surface.




Backrest handle

Small handle on the side of the backrest, convenient for folding and unfolding the chair, less effort.

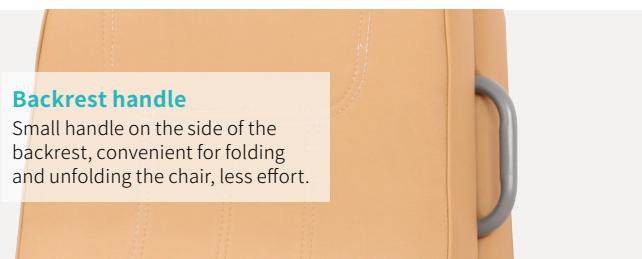

Bed foot

Bed legs are equipped with plastic cover to prevent possible noise caused by friction with the floor surface; when folded the backrest is securely fixed, to enhance patient's safety.

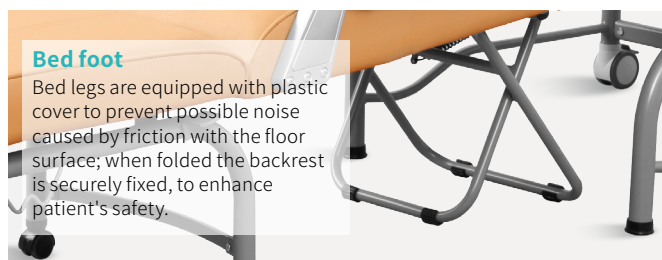

Function

Telescopic springs inside, material 65Mn, self-positioning design, two easy switchable modes.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd
Sala Thammassop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited
DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya Naberezhnaya,
Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue