



SKE942

RECLINER CHAIRS FOR
ELDERLY



SAIKANG / Medical

Healthier
your life

SKE942

FEATURES AND HIGHLIGHTS

Technical parameters

▪ External size	1205×668×1180mm
▪ Seat size	530×500×520mm
▪ Tray size	650×390mm
▪ Back-rest adjustment	37°-61°
▪ Leg-rest adjustment	13°-36°
▪ Safe working load	113kg



Design

Manual recliner with four sections (back, hip, leg, foot), ergonomic design, ideal for dialysis, blood collection, respiratory and cardiac care, and adapts to patient needs.



Mechanical Function

Manual adjustment of the backrest, three modes: sitting, semi-lying, lying; footrest is completely hidden when the chair in the sitting mode, to prevent bumps and falls.



Surface material

Using shaped foam sponge, PU leather cover, handmade finishing. Option color: red/ green/ blue.




Detachable dining table

The dining table board adopts the plug-in installation, with 5 positions, it can be adjusted through the buckle, suitable for different patients, and can prevent demented patients from leaving the seat at will and falling.



Caster

5-inch caster with brake, easy to drive, the loading capacity up to 113kg.

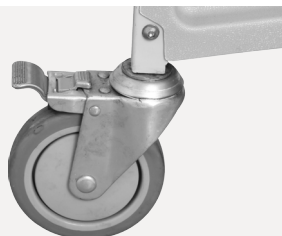



CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd
Sala Thammassop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited
DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya Naberezhnaya,
Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue