



A106-2

ELECTRIC-HYDRAULIC OPERATING TABLE

Registration No. SIF-B



SAIKANG / Medical

Healthier
your life

A106-2

FEATURES AND HIGHLIGHTS



Technical parameters

External Length:	2080mm
External Width:	550mm
Table Height Min:	500mm
Table Height Max:	800mm
Reverse Trendelenburg:	0°-25°
Trendelenburg:	0°-25°
Lateral Tilt Left/Right:	0°-20°/0°-20°
Head Plate Up/Down:	0°-45°/0°-90°
Back Plate Up/Down:	0°-80°/0°-20°
Foot Plate Up/Down:	0°-15°/0°-90°

Foot Plate Horizontal Positioning:	0°-90°
Patient Weight Capacity:	250kg
X-rays translucent:	YES
C-arm compatible:	NO
Multi-functional orthopedic frame:	YES

Technical configuration

<input checked="" type="checkbox"/> Anaesthesia Screen	1pc
<input checked="" type="checkbox"/> Hand Remote	1pc
<input checked="" type="checkbox"/> Kidney Bridge Handle	1pc
<input checked="" type="checkbox"/> Arm Support	1pc
<input checked="" type="checkbox"/> Foot Support	1pc
<input checked="" type="checkbox"/> Mattress	1set
<input checked="" type="checkbox"/> Shoulder Support	1pair
<input checked="" type="checkbox"/> Long Fixing Clamp	1pair
<input checked="" type="checkbox"/> Fixing Clamp	7pcs
<input checked="" type="checkbox"/> Power Line	1pc

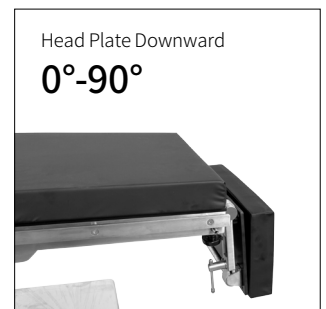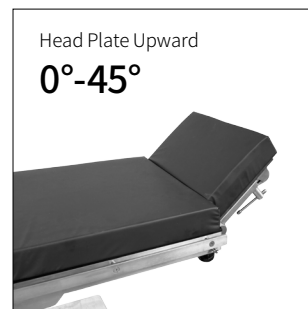

Material

The base, column and side hoder are all made of high quality 304 stainless steel, which is anti-corrosive and can be easily cleaned.



Powder Supply

Internal battery maintaining power supply for 50-80 operations and with low battery warning function.




Remote Panel

All function controlled by hand controller in the column, and can be locked.



Hydraulic System

Brand hydraulic system (German Hoerbiger) make table with low noise, high safety and reliability.

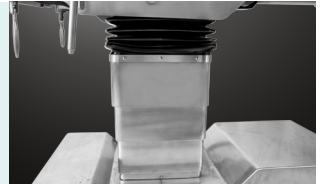

Control System

Foot treadle brake device which is easily for moving operating table.



Super Low Position

Super low position (500mm) , suitable for operation on head.



Kidney Bridge

Built-In Kidney Bridge, More Convenient And Nice Appearance.



Screen



Armrest



Knee Curtch



Shoulder Support



Body Support





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc

W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd

Sala Thammassop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited

DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'

454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC

Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue