



A2000

ELECTRIC OPERATING TABLE

Registration code: SKL-A



SAIKANG / Medical

Healthier
your life

A2000

FEATURES AND HIGHLIGHTS



Technical parameters


External size (LxWxH):	2100x480x(680-930)mm
Reverse trendelenburg:	0-25°
Trendelenburg:	0-25°
Lateral tilt left/right:	0-15°
Head plate up/down:	45°/90°
Back plate up/down:	70°/15°
Leg plate up/down:	45°/90°
Leg plate outwards:	0-90°
Kidney bridge elevation:	0-100mm
Safe working load:	150kg
X-rays translucent:	NO
C-arm compatible:	NO
Multi-functional orthopedic frame:	YES

Technical configuration

<input checked="" type="checkbox"/> Bed frame	1 pc
<input checked="" type="checkbox"/> Mattress	1 set
<input checked="" type="checkbox"/> Anesthesia screen frame	1 pc
<input checked="" type="checkbox"/> Shoulder holder	2 pcs
<input checked="" type="checkbox"/> Arm holder	2 pcs
<input checked="" type="checkbox"/> Leg support	2 pcs
<input checked="" type="checkbox"/> Foot support	2 pcs
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Sliding block	10 pcs
<input checked="" type="checkbox"/> Spanner	4 pcs


Material

The base material for the table is high thickness 304 stainless steel, seamless welding by robot, with modern appearance and more stable.

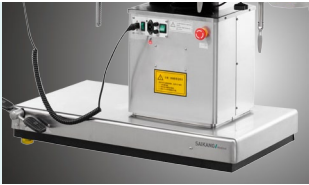


Mattress

Removable memory mattress, good rebound resilience, anti-static electricity and anti-microbial.

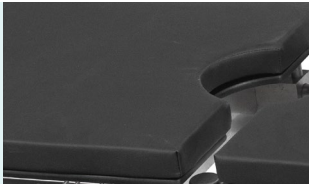

Brake

Pedal brake, easy to move and lock.

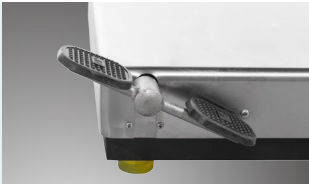

Safety

The table is provided with an emergency stop switch on the base column, help to stop the operating table running in case of emergency situation, ensure the safety of operation.



Leg plate

Leg plate with imported gas spring for easy control.




Kidney bridge

Manual kidney bridge, more convenience.



Battery could support

Battery could support about 50-80 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.

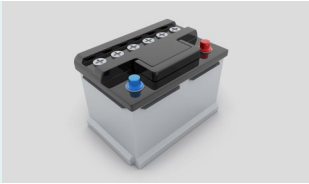

Armrest



Anesthesia Frame




Leg Crutches



Shoulder Support



Body Supporter

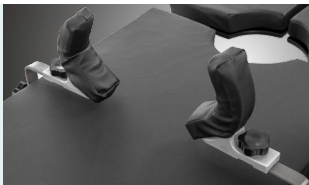

Professional functions

Lowest position of table top

680mm



Highest position of table top

930mm



Back plate up

0°-70°



Back plate down

0°-15°



Reversed trendelenburg

0°-25°



Trendelenburg

0°-25°



Right lateral tilt

0°-15°



Left lateral tilt

0°-15°



Head plate upward

0°-45°



Head plate downward

0°-90°



Leg section evagination

0°-90°



Leg plate up/down

-45°-90°



Kidney bridge elevation

0-100mm





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc

W3400 Corrientes, Argentina
 TEL: +55 (11) 2758-0438
 Rua Universal, Hollywood,
 São Bernardo do Campo,
 Brazil
 TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
 upper hill, nairobi, Kenya
 TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd

Sala Thammasop Subdistrict
 Thawi
 Watthana District
 Bangkok Thailand
 TEL: +66 642741046

UK Branch

Uk Healthward International Limited

DEPT 111,
 196 High Road Wood Green,
 London, N22 8HH
 United Kingdom
 TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'

454004, Russian Federation,
 Chelyabinsk,
 Universitetskaya
 Naberezhnaya, Russia
 TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC

Musaffah
 P.O Box 53200, Abu Dhabi,
 UAE
 Tel: +971 585710652

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
 STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue