



# A301

## ELECTRIC OPERATING TABLE

Registration code: SKL-A



**SAIKANG** / Medical

Healthier  
your life

# A301

## FEATURES AND HIGHLIGHTS



### Technical parameters


External size (LxWxH):	2150x555x(750-1000)mm
Reverse trendelenburg:	0-25°
Trendelenburg:	0-25°
Lateral tilt left/right:	0-15°
Head plate up/down:	30°/90°
Back plate up/down:	70°/15°
Leg plate down:	90°
Leg plate outwards:	90°
Kidney bridge elevation:	0-110mm
Horizontal sliding:	0-350mm
Safe working load:	275kg
X-rays translucent:	YES
C-arm compatible:	YES
Multi-functional orthopedic frame:	YES

### Technical configuration

<input checked="" type="checkbox"/> Bed frame	1pc
<input checked="" type="checkbox"/> Mattress	1set
<input checked="" type="checkbox"/> Anesthesia screen frame	1pc
<input checked="" type="checkbox"/> Shoulder holder	2pcs
<input checked="" type="checkbox"/> Arm holder	2pcs
<input checked="" type="checkbox"/> Leg support	2pcs
<input checked="" type="checkbox"/> Foot support	2pcs
<input checked="" type="checkbox"/> Hand controller	1pc
<input checked="" type="checkbox"/> Sliding block	8pcs
<input checked="" type="checkbox"/> Spanner	5pcs

**Material**

The base material for the table is high thickness 304 stainless steel, seamless welding by robot, with modern appearance and more stable.

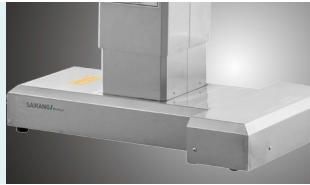


**Mattress**

Polyurethane foam mattress, waterproof and seamless, avoid liquid in. The connection between table and mattress is buckle joint, easy to disassemble and install, easy to clean.



**Head plate**

Head plate is with imported gas spring for easy control.

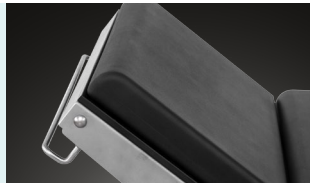

**Kidney bridge**

Manual kidney bridge, operate by head side rocker, more convenience.



**Horizontal sliding**

Extra horizontal sliding, distance 350mm high strengthen phenolic resin table surface, C-arm and X-ray compatible.



**Dual control system**

All functions controlled by both hand controller or control panel on table base column. With action LCD display on hand controller, easy to operate. With control panel on table base column, in case of emergency hand controller failure.



**Leg plate**

Leg plate with imported gas spring for easy control. Disassembled type leg plate.



**Brake**

Pedal brake, easy to move and lock. Pedal down to release the inner castors down. Pedal up to retract the castors, which enables the table locked on the floor.



**Battery could support**

Battery could support about 80-120 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.



**Armrest**




**Anesthesia Frame**



**Leg Crutches**



**Shoulder Support**

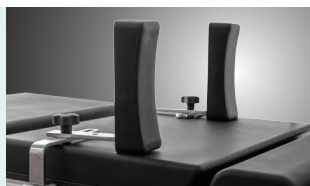

**Body Supporter**



Professional functions

Lowest position of table top  
**750mm**



Highest position of table top  
**1000mm**



Back plate up/down  
**-15°-70°**



Kidney bridge elevation  
**0-110mm**



Reversed trendelenburg  
**0°-25°**



Trendelenburg  
**0°-25°**



Right lateral tilt  
**0°-15°**



Left lateral tilt  
**0°-15°**



Head plate upward  
**0°-30°**



Head plate downward  
**0°-90°**



Leg section evagination  
**0°-90°**



Leg plate down  
**0°-90°**





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### South America Branch

#### Saikang Medical Technology, Inc

W3400 Corrientes, Argentina  
TEL: +55 (11) 2758-0438  
Rua Universal, Hollywood,  
São Bernardo do Campo,  
Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

#### Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

#### Saikang Medical Technology Ltd

Sala Thammasop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

#### Uk Healthward International Limited

DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 202111106

### Russia Branch

#### LLC 'SAIKANG'

454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya  
Naberezhnaya, Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

#### Saikang Medical Supplies (Mid East) LLC

Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue