



A302

ELECTRIC
OPERATING TABLE

Registration code: SKL-A

SAIKANG / Medical

Healthier
your life

A302

FEATURES AND HIGHLIGHTS



Technical parameters

External size (LxWxH):	2150x555x(700-1000)mm
Reverse trendelenburg:	0-25°
Trendelenburg:	0-25°
Lateral tilt left/right:	0-25°
Head plate up/down:	45°/90°
Head plate Height adjustment:	70mm
Back plate up/down:	70°/22°
Leg plate up/down:	15°/90°
Leg plate outwards:	90°
Kidney bridge elevation:	0-110mm
Horizontal sliding:	0-350mm
Safe working load:	275kg
X-rays translucent:	YES
C-arm compatible:	YES
Multi-functional orthopedic frame:	YES

Technical configuration

<input checked="" type="checkbox"/> Bed frame	1pc
<input checked="" type="checkbox"/> Mattress	1set
<input checked="" type="checkbox"/> Anesthesia screen frame	1pc
<input checked="" type="checkbox"/> Shoulder holder	2pcs
<input checked="" type="checkbox"/> Arm holder	2pcs
<input checked="" type="checkbox"/> Leg support	2pcs
<input checked="" type="checkbox"/> Foot support	2pcs
<input checked="" type="checkbox"/> Hand controller	1pc
<input checked="" type="checkbox"/> Sliding block	8pcs
<input checked="" type="checkbox"/> Spanner	5pcs

Material

The base material for the table is high thickness 304 stainless steel, seamless welding by robot, with modern appearance and more stable.




Mattress

Polyurethane foam mattress, waterproof and seamless, avoid liquid in. The connection between table and mattress is buckle joint, easy to disassemble and install, easy to clean.



Head plate

Head plate is with imported gas spring for easy control.

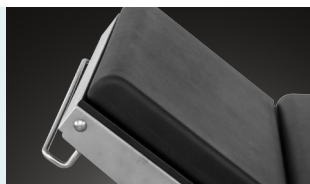

Kidney bridge

Manual kidney bridge, operate by head side rocker, more convenience.



Horizontal sliding

Extra horizontal sliding, distance 350mm high strengthen phenolic resin table surface, C-arm and X-ray compatible.



Dual control system

All functions controlled by both hand controller or control panel on table base column. With action LCD display on hand controller, easy to operate. With control panel on table base column, in case of emergency hand controller failure.



Leg plate

The up and down folding actions of the leg plate are all completed by electric. The leg plate is removable, effortless and smooth.



Castors

HAION brand castors, each load bearing is 250kg, smooth and noiseless movement, with longer service life. Pedal brake, easy to move and lock. Unique switch design.




Battery could support


Battery could support about 80-120 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.



Armrest

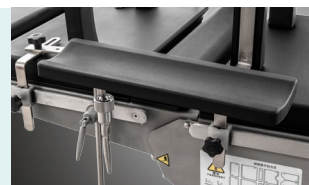


Anesthesia Frame

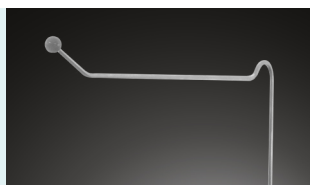

Leg Crutches



Shoulder Support

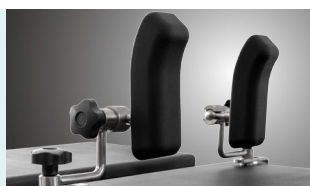

Body Supporter



Professional functions

<p>Lowest position of table top 700mm</p> 	<p>Highest position of table top 1000mm</p> 	<p>Back plate up 0°-70°</p> 	<p>Back plate down 0°-22°</p> 
<p>Reversed trendelenburg 0°-25°</p> 	<p>Trendelenburg 0°-25°</p> 	<p>Right lateral tilt 0°-25°</p> 	<p>Left lateral tilt 0°-25°</p> 
<p>Head plate upward 0°-45°</p> 	<p>Head plate downward 0°-90°</p> 	<p>Leg section evagination 0°-90°</p> 	<p>Leg plate up/down -15°-90°</p> 
<p>Head plate Height adjustment 0-70mm</p> 	<p>Kidney bridge elevation 0-110mm</p> 	<p>Horizontal sliding 0-350mm</p> 	
			



CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc

W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd

Sala Thammasop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited

DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'

454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC

Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue