



# SKB041-3

## PATIENT TRANSPORTATION TROLLEY

Registration code: SKM-F




**SAIKANG** Medical

Healthier  
your life

# SKB041-3

## FEATURES AND HIGHLIGHTS

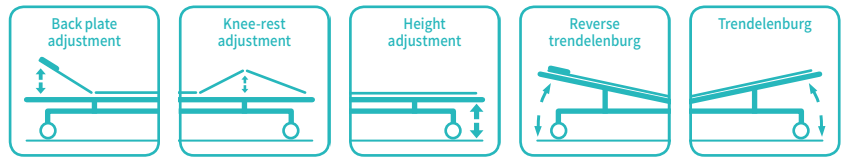

### Technical parameters

#### X-ray Backrest

▪ External size (LxW):	2120x840 mm
▪ Height from floor to bed platform:	550-860 mm
▪ Height from floor to lowered guardrail:	590-900 mm

#### X-ray Platform

▪ External size (LxW):	2120x840 mm
▪ Height from floor to bed platform:	600-910 mm
▪ Height from floor to lowered guardrail:	640-950 mm

▪ Back plate adjustment:	0-85° (±5°)
▪ Knee-rest adjustment:	0-30° (±5°)
▪ Reverse Trendelenburg:	0-15°
▪ Trendelenburg:	0-15°

### Technical configuration

<input checked="" type="checkbox"/> Usa hydraulic pump	2 pcs
<input checked="" type="checkbox"/> Crank system	1 set
<input checked="" type="checkbox"/> Back rest air spring	2 pcs
<input checked="" type="checkbox"/> 8" German castors,central locking	4 pcs
<input checked="" type="checkbox"/> Emergency push handle	1 set
<input checked="" type="checkbox"/> Directional wheel	1 pc
<input checked="" type="checkbox"/> Folding side rail	1 set
<input checked="" type="checkbox"/> Φ74mm wall bumper wheel	4 pcs
<input checked="" type="checkbox"/> Iv pole previsions	4 pcs
<input checked="" type="checkbox"/> Mattress	1 pc
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> Central locking break pedal	2 sets
<input checked="" type="checkbox"/> Directional wheel release pedal	2 pcs
<input checked="" type="checkbox"/> X-ray translucent bed platform	1 set



**Metal material**

Baosteel from Fortune Global 500  
Panasonic Robotic ensure 360° full  
smooth wedding.

**Plastic**

Environmental plastic can be  
sterilized in 100°C, tensile upto  
30MPa.

**Steels | Painting**

11 Process epoxy painting, ASTM  
testing anti-bacterial, paint thickness  
0.12 mm, brightness 60°, paint can  
resist 50 kg impact.

**Brake**

Break spring keep qualified after  
10,000times uninteruped break.

**Hydraulic pumps**


Height adjustment, trendelenburg  
and revised tendelenburg through  
two USA hydraulic pumps.

**Center wheel**

Can be engaged to facilitate direction  
control, enabling users to navigate  
the stretcher easily in hectic and  
dynamic environments.

**Rising**

Hydraulic treadle for rise.

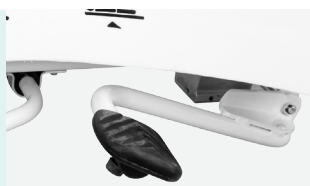**Falling**

① Leg section lows down to make  
head higher.  
② Back section lows down to make  
leg section higher.  
③ Step on the middle of the pedal,  
make the bed base fall down steadily.

**Castors**

Four double side casters of 8 inch  
diameter, simultaneously lockable by  
the step pedal at both sides.  
TPR tire no worn out after running  
30KM, save anti-winding hard shell,  
united forming without bolts.

**Mattress**

Waterproof PU mattress, convenient  
for transfer patient.

**Side rail**

Foldable aluminium alloy guardrail

**Bumper wheel**


Bumper wheel detachable, max.  
protect the bed and patient during  
transportation

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic  
hook.

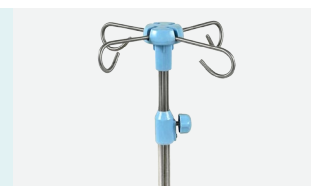





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### South America Branch

#### Saikang Medical Technology, Inc

W3400 Corrientes, Argentina  
TEL: +55 (11) 2758-0438  
Rua Universal, Hollywood,  
São Bernardo do Campo,  
Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

#### Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

#### Saikang Medical Technology Ltd

Sala Thammasop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

#### Uk Healthward International Limited

DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 202111106

### Russia Branch

#### LLC 'SAIKANG'

454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya  
Naberezhnaya, Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

#### Saikang Medical Supplies (Mid East) LLC

Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue